

### BZ-G03 Green™ Portable oxygen generator



Features


- \* Continuous oxygen supply for 24 hours/day
- \* High performance, stability and portability
- \* Low noise up to 46 dB for 3L oxygen concentrator
- \* Automatic time accumulation function (max 15,000 hours) and automatic alarm for low oxygen concentration
- \* High performance, energy-saving advanced Lithium molecular sieve from Switzerland
- \* User friendly luggage pull handle design
- \* Can be used indoors, at nursing homes, in hospitals, office rooms and other places

## Imported molecular sieve



## PSA(Pressure swing adsorption)


The filtered air will go into the compressor and it will compress the air into the molecular sieves, one of the molecular sieve produces oxygen with PSA and another starts low pressure desorption. The oxygen will then store in oxygen tank before gives to patient.

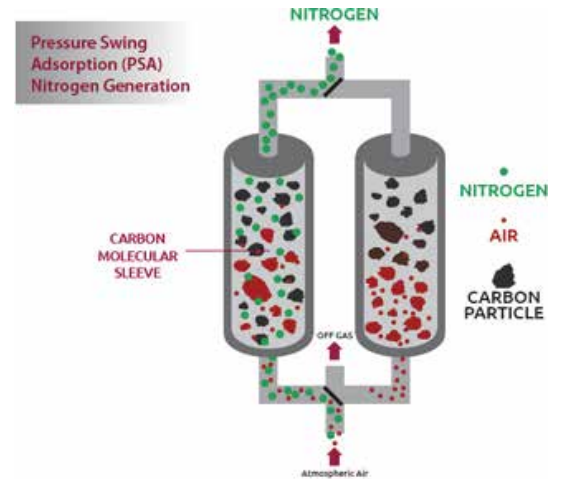

## Nebulization Function

- One - machine, Dual - use
- Nebulizer makes the mist or liquid medicine into bronchus and pulmonary alveoli, treatment for most of the respiratory disease.


## Specification

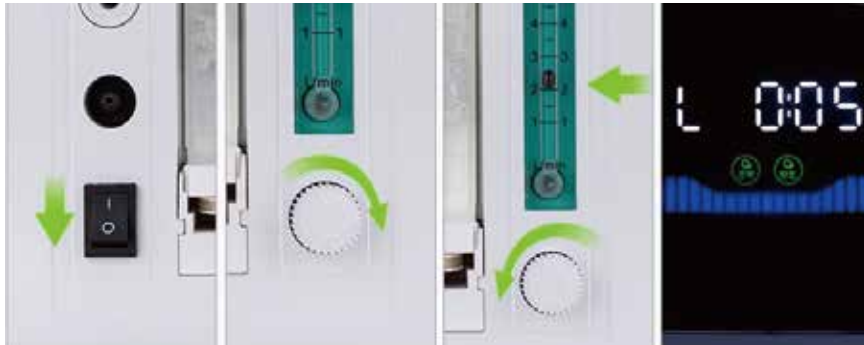

1. Turn on the power switch
2. Rotate the flow knob to the end of right.
3. Rotate left stop the floater between 2 and 3.
4. Oxygen supply begins after 5 minutes

Model Num.	BZ-G03
Type	Home/Medical Oxygen Concentrator
Flow	1-5L/min
Voltage	220±22V, 50±1 HZ
Valve Plate	Switzerland
Rated Power(VA)	≤230
Oxygen Concentration	>90%
Noise	≤46dB
Operating Hours	24 hours/day
Type of Molecular Sieve	Lithium (Switzerland)
Output Pressure	0.04-0.05Mpa
Display	LCD
Size	300*200*450mm
Weight	8.8kg

High-purity green oxygen  
— Multiple filtration system

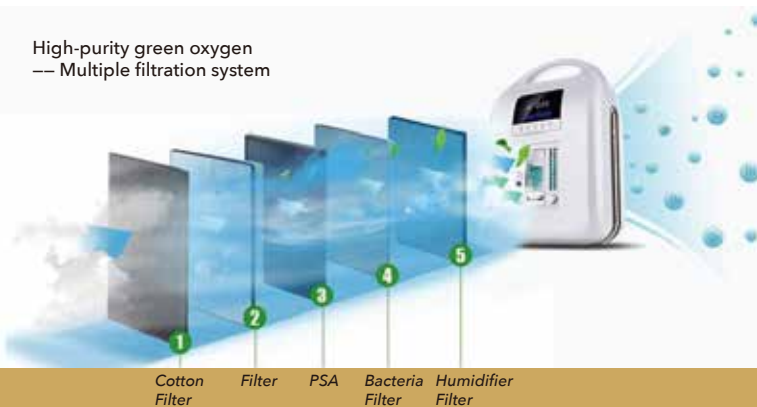

Cotton Filter    Filter    PSA    Bacteria Filter    Humidifier Filter



8.8 kg  
Integrated portable handle design

### Accessories

